


From start to finish, a PMD is convenient and easy to use.



- 1 Medication Cups** Comfortable to hold with secure lids to prevent spills. Each reusable cup can hold 20-25 pills and supplements. As many as six cups/doses can be delivered per day.
- 2 Dispensing Button** Designed for seniors' special needs. Easy to see and feel, it is conveniently oversized and requires only a gentle push for one-button dispensing.
- 3 Medication Reminders** Each dose is announced by a loud clear voice and tone (adjustable volume), a visual message, and a flashing red light. Special instructions such as "Take with food" can also be announced.
- 4 Keypad** Designed for easy setup, it's located in the locking compartment to prevent accidental changes or tampering.
- 5 Battery Backup** Provides alternate power for up to 18 hours.
- 6 Locking Storage Compartment** It can be pre-loaded with a maximum of 60 medication cups. The locking compartment keeps medications safe and in order.

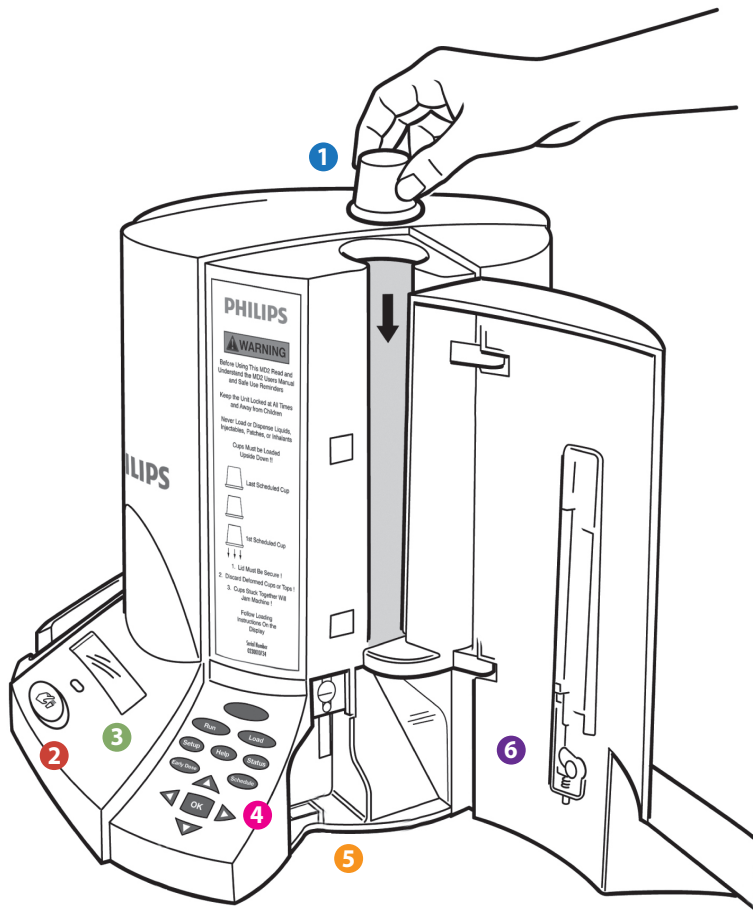

1

To get you started, H&HC will help you set up the medication schedule right over the phone.



2

Filling and refilling the medication cups is easy. Just use one cup for morning pills, one cup for afternoon pills, etc. An organizing tray helps you load cups in the right order.



3

The PMD announces when it's time to take each dose. Then, just pushing a button delivers the right pills.



4

If the medications aren't dispensed after multiple reminders, the PMD system automatically alerts family members or friends by telephone.



5

If medications run low, the PMD will alert you by phone that it's time to re-load.